







Leibniz-Institut für Troposphärenforschung Permoserstraße 15 04318 Leipzig

Intercomparison of Condensation Particle Counter

Project No.: CPC-2019-3-10

Principal Investigator: Kay Weinhold

Home Institution: TROPOS

Participant:

Candidate: CPC Ref1

Counter (SN): TSI CPC Model 3772 SN3772141701

Location of the quality assurance: TROPOS Leipzig, lab 130

Comparison period: July 10, 2019

Last Intercomparison (with Project No.):

TROPOS Reference Instrument: Electrometer: TSI model 3068B

#70838596, Last calibration in September 2018

Additional Equipment: Bubble flow meter 'Gilibrator', Gilian (Sensidyne)

#1711008-S, Last calibration in January 2018

Summary of Intercomparison

Status:

The candidate passed the quality standards of ACTRIS and GAW. The candidate reached 100% efficiency at 40 nm. The Dp50 is at 8.85 nm. The CPC efficiency curve corresponds to the standard of ACTRIS and GAW.

Page 1 / 4








for Aerosol Physics





eibniz Institute for Tropospheric Research

Leibniz-Institut für Troposphärenforschung Permoserstraße 15 04318 Leipzig

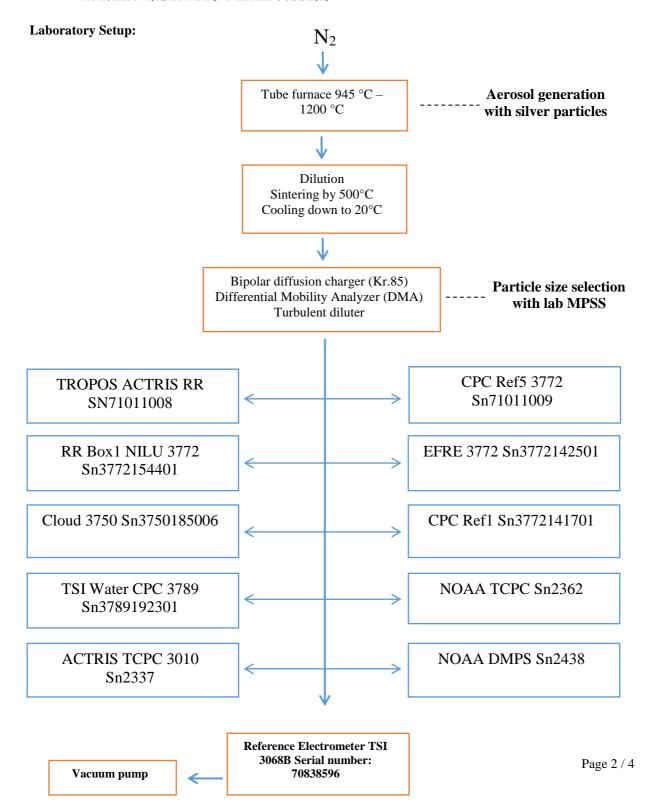

Leibniz-Institut für Troposphärenforschung e.V.
Telefon: +49 341 2717-7060

Telefax: +49 341 2717-99-7060 info@tropos.de http://www.tropos.de Commerzbank Leipzig KTO 102 14 50 BLZ 860 400 00 IBAN: DE77 8604 0000 0102 1450 00

SWIFT CODE: COBADEFF 860

Leibniz-Gemeinschaft









Leibniz Institute for Tropospheric Research

Leibniz-Institut für Troposphärenforschung Permoserstraße 15 04318 Leipzig

Date of arrival of instrument in calibration lab: July 10, 2019

Instrument:

Condensation Particle Counter Model and serial number of instrument: CPC 3772 S/N 3772141701

for Aerosol Physics

Result of physical inspection: no damages

Result of functional test: functional test successful

Internal parameters of instrument nominal flow rate 1.0 l/min

Model and identification number of

aerosol electrometer: TSI Electrometer Model 3068, S/N 70838596

Electrometer calibration certificate: September 05, 2018, calibrated at PTB

Braunschweig

Corrections of electrometer, for instance,

differing flow rate:

Within tolerance range (+/-2%); reference: 4.0

l/min, measured: 3.970 l/min

Software for recording: LabView 2010; National Instruments; Program

"LabCount.vi"

Date of calibration: July 10, 2019

Lab temperature and pressure: 23°C, 983.82mbar Measured aerosol flow rate of CPC: 1.017 l/min

Uncertainty in measured flow rate: 3%

Flowmeter used: Gilian Gilibrator V; S/N 1711008-S,

January, 2018

Particles and gases used for calibration: silver particles and nitrogen Method of particle generation: tube furnace generator **Zero measurement of instrument:** 0 particles/cm³ in 5 minutes


Results (using pulse output):

Particle size (nm)	40	30	20	15	12
Number concentration (cm-3)	1221	1327	1366	1492	1053
Counting efficiency η	1.01	1.01	0.99	0.92	0.82
Particle size (nm)	10	09	08	07	06
Number concentration (cm-3)	859	724	604	325	11
Counting efficiency η	0.66	0.52	0.35	0.15	0.01
Particle size (nm)	40				
Number concentration (cm-3)	1194				
Counting efficiency η	1.01				

SWIFT CODE: COBADEFF 860





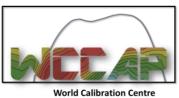

for Aerosol Physics





Leibniz Institute for Tropospheric Research

Leibniz-Institut für Troposphärenforschung Permoserstraße 15 04318 Leipzig

Special Information regarding to the Candidate:

Was it necessary to:	yes/no	information
do a second run	no	-
clean the optics	no	-
clean the nozzle	no	-
clean the saturator	no	-
change the wick	no	-
change the laser	no	-
change internal settings	no	-

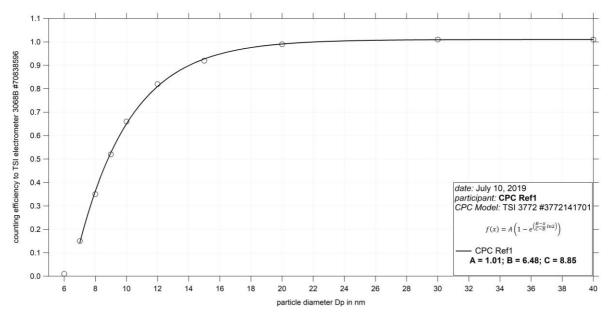


Fig. 1: Counting efficiency for CPC Ref1 3772 SN3772141701 against aerosol electrometer 3068 SN 70838596; silver particles between 6 and 40 nm were used for calibration; the calculated Dp50 is 8.85 nm.

Status information:

244445 1111011114110111					
Status	T SAT	T CON	T OPT	T CAB	P AMB
from display	39.0	22	40.0	29.6	99.6
Status	P OR	P NO	Laser	LV	flow
from display	82.6	2.8	50	full	1.017

Date of issue: July 10, 2019

Reviewed: TROPOS / Kay Weinhold Page 4 / 4



SWIFT CODE: COBADEFF 860